
HYDRA FRESH (WFA340)

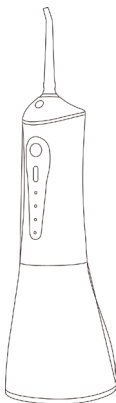ES. Instrucciones de uso
EN. Instructions for use
FR. Mode d'emploi
PT. Instruções de uso
IT. Istruzioni per l'uso
CA. Instruccions d'ús
DE. Gebrauchsanleitung
NL. Gebruiksaanwijzing
RO. Instrucțiunile de utilizare
PL. Instrukcja obsługi

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ORAL IRRIGATOR WFA340

Dear customer,

Many thanks for choosing to purchase a TAURUS brand product.

Thanks to its technology, design and operation and the fact that it exceeds the strictest quality standards, a fully satisfactory use and long product life can be assured.

SAFETY ADVICE AND WARNINGS

Read these instructions carefully before switching on the appliance and keep them for future reference. Failure to follow and observe these instructions could lead to an accident.

Oral irrigators can be used by children aged from 8 years

and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.

This oral irrigator shall be used only with cold or warm water (below 40°C).

Children shall not play with the appliance.

This appliance is for household use only, not professional or industrial use.

For the purposes of recharging the battery, only use a Safety Extra Low Voltage power supply (SELV).

This appliance contains batteries that are non-replaceable. When the battery is

at end of life, the appliance shall be properly disposed of.

Ensure that the voltage indicated on the rating label matches the mains voltage before plugging in the appliance.

The device battery can only be replaced by specialized personnel from an authorized technical support service.

Do not use the device when the USB cable or connector is damaged.

Do not wrap the USB cable around the appliance.

If the appliance enclosure breaks, immediately disconnect the appliance from the mains to prevent the possibility of an electric shock.

Do not touch the plug with wet hands.

DESCRIPTION

- A. Nozzle release button
- B. ON/OFF button
- C. Mode button
- D. Main body
- E. Pilot light
- F. Water tank

- G. Water inlet lid
- H. Charging port
- I. Nozzles
 - K1. Standard nozzle (x1)
 - K2. Dental brush nozzle (x1)
 - K3. Tongue cleaner nozzle (x1)
 - K4. Periodontal nozzle (x1)
- J. USB cable
- K. Storage bag

If the model of your appliance does not have the accessories described above, they can also be bought separately from the Technical Assistance Service.

USE AND CARE:

- Do not use the appliance if the parts or accessories are not properly fitted.
- Do not use the appliance if the accessories attached to it are defective. Replace them immediately.
- Do not use the appliance if it is empty.
- Do not use the appliance if the ON/OFF button does not work.
- Do not use the appliance if it is tipped up and

do not turn it over.

- Do not turn the appliance over while it is in use or connected to the charger.
- If you use the appliance in a bathroom or similar place, unplug the appliance from the mains when it is not in use, even if it is only for a short while, as proximity to water involves a risk, even if the appliance is disconnected.
- Disconnect the USB cable from the mains when not in use and before undertaking any cleaning task.
- Keep the appliance in a dry, dust-free place, out of direct sunlight.
- Keep the appliance in good condition. Check that the moving parts are not misaligned or jammed and make sure there are no broken parts or anomalies that may prevent the appliance from operating correctly.
- Use the appliance and its accessories and tools in accordance with these instructions, taking into account the working conditions and the work to be performed. Use the appliance for operations different from those intended could result in a hazardous situation.
- Never leave the appliance connected and unattended if it is not in use. This saves energy

and prolongs the life of the appliance.

- Do not use the appliance for more than 8 minutes at a time. When using exceeds 8 minutes cycles allow a resting period of at least 2 minutes between each cycle. Under no circumstances should the appliance be operated for more time than necessary.
- Do not use the appliance with pets or animals.

CHARGING THE BATTERY

- Caution: A new battery is not fully charged and must therefore be charged before the appliance is used for the first time.
- Charge the battery at an atmospheric temperature of between 4 and 40°C.
- Connect the USB cable to the mains or adapter.
- A pilot light will come on to indicate that the charger is ready to begin charging the battery.
- Approximately 6 hours are needed for a complete charge of the battery. Never attempt to over charge the battery (although it is recommended that the first charge of the battery be of 6 hours in duration).
- Insert the USB cable to the charging port. The charging port is located at the back part of the

appliance.

- The pilot light will come on to indicate that the charger is charging the battery.
- When the battery is charged disconnect the USB from the mains supply.
- NB: Do not use the appliance while it is connected to the charger.

HOW TO CHECK THE BATTERY'S CHARGE LEVEL:

- Charging: The white charging indicator light flashes slowly.
- Low battery: The white charging indicator light flashes.
- Fully charged: The white charging indicator light stays on.

INSTRUCTIONS FOR USE BEFORE USE:

- Make sure that all the product's packaging has been removed.
- Please read carefully the "Safety advice and warnings" section before first use
- Before using this product for the first time, it is important to make a full charge of battery.

FILLING WITH WATER:

- You must fill the water tank before turning on the appliance.
- Open the water inlet lid.
- Fill the water tank with cold water.
- Close the water inlet lid.

INSERTING AND REMOVING NOZZLES

- Insert the nozzle into the appliance and press down until you hear a “click”.
- To remove the nozzle from the appliance, press the nozzle release button and pull out the nozzle.

SPRAY NOZZLES

STANDARD NOZZLE:

- Start rinsing from the molars towards the incisors.
- Let the spray head lightly shift along the gums, and pause briefly between teeth.
- Repeat with each side of the teeth.

DENTAL BRUSH NOZZLE:

- Place the brush head directly against the teeth so the bristles touch the enamel and gumline.

- Manually brush each tooth with a light scrubbing motion to dislodge plaque while the water washes it away.
- Guide the nozzle slowly from the back teeth to the front, cleaning both the surface and the spaces between.

TONGUE CLEANER NOZZLE:

- Lean low over your sink, place the tongue cleaner on your tongue and start the appliance on the lowest pressure setting, increase as comfortable.
- Place the tongue cleaner in the center of your tongue about half way back and pull forward with light pressure.

PERIODONTAL NOZZLE

- The periodontal nozzle delivers water or anti-bacterial solutions deep into periodontal pockets.
- Select the soft mode.
- Lean low over your sink and place the nozzle against a tooth at a -45degree angle and gently place the tip under the gum line inserting the nozzle between teeth.
- Continue tracing along the gum line until you finish cleaning.

USE:

- Select the nozzle you are willing to use and place it until you hear a “click”.
- Guide the appliance to direct the flow of water in the desired direction.
- Turn the appliance on, by using the ON/OFF button.
- Select the desired mode using the mode button.

ONCE YOU HAVE FINISHED USING THE APPLIANCE:

- Turn the appliance off pressing the ON/OFF button.
- Remove the water from the water tank.
- Clean the appliance.

SAFETY BATTERY PROTECTOR:

- The appliance has a lithium battery protection It will automatically down off in case of overcurrent or short circuit, and can be used again after connecting to power.

CLEANING

- Clean the equipment with a damp cloth with a

few drops of washing-up liquid and then dry.

- Do not use solvents, or products with an acid or base pH such as bleach, or abrasive products, for cleaning the appliance.
- Never submerge the appliance in water or any other liquid or place it under running water.
- Never submerge the connector in water or any other liquid or place it under running water.

HOW TO REMOVE THE BATTERY FROM THE INTERIOR OF THE DEVICE

- This appliance contains battery that is only replaceable by skilled persons.
- In case you need to replace the battery, take the appliance to an authorised technical support service.
- Do not try to dismantle or repair without assistance, as this may be dangerous.

SERVICE:

Any misuse or failure to follow the instructions for use renders the guarantee and the manufacturer's liability null and void.

ANOMALIES AND REPAIR

Take the appliance to an authorised technical support service if problems arise. Do not try to dismantle or repair without assistance, as this may be dangerous.

FOR EU PRODUCT VERSIONS AND/OR IN THE CASE THAT IT IS REQUESTED IN YOUR COUNTRY:

ECOLOGY AND RECYCLABILITY OF THE PRODUCT

The materials of which the packaging of this appliance consists are included in a collection, classification and recycling system. Should you wish to dispose of them, use the appropriate public containers for each type of material.

The product does not contain concentrations of substances that could be considered harmful to the environment.

This symbol means that in case you wish to dispose of the product once its working life has ended, take it to an  authorised waste agent for the selective collection of waste electrical and electronic equipment (WEEE).

This symbol means that product may include a battery or batteries; user must

to remove them before disposing of the product. Remember that batteries must be disposed of in duly authorised containers. Do not dispose of them in fire.

This appliance complies with Directive 35/2014/EU on Low Voltage, Directive 30/2014/EU on Electromagnetic Compatibility, Directive 65/2011/EU on the restrictions of the use of certain hazardous substances in electrical and electronic equipment and Directive 125/2009/EC on the ecodesign requirements for energy-related products.

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ENGLISH

WARRANTY AND TECHNICAL ASSISTANCE

- This product enjoys the recognition and protection of the legal guarantee in accordance with current legislation. To enforce your rights

or interests you must go to any of our official technical assistance services.

- You can find the closest one by accessing the following web link: <http://taurus-home.com/>
- You can also request related information by contacting us.
- You can download this instruction manual and its updates at <http://taurus-home.com/>

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