
HYDRA SMILE

- ES. Instrucciones de uso
- EN. Instructions for use
- FR. Mode d'emploi
- PT. Instruções de uso
- IT. Istruzioni per l'uso
- CA. Instruccions d'ús
- DE. Gebrauchsanleitun
- NL. Gebruiksaanwijzing
- RO. Instrucțiuni de utilizare
- PL. Instrukcja obsługi
- BG. Инструкция за употреба
- EL. Οδηγίες χρήσης
- RU. Инструкция по применению
- DA. Brugsanvisning
- NO. Instruksjoner for bruk
- SV. Användningsinstruktioner
- FI. Käyttöohjeet
- TR. Kullanım için talimatlar
- HE. הוראות לשימוש
- AR. تعليمات الاستخدام

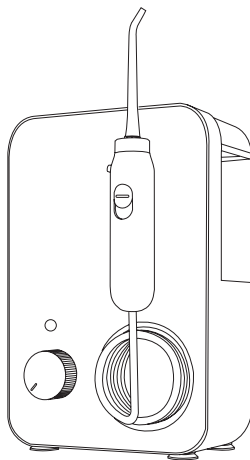taurus

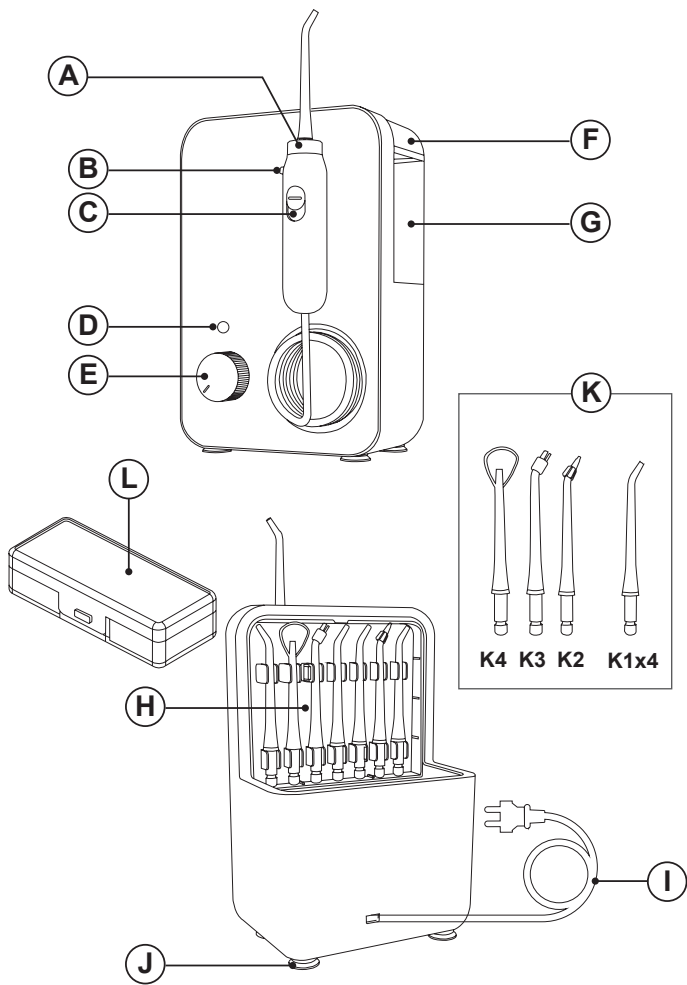
Fig. 1


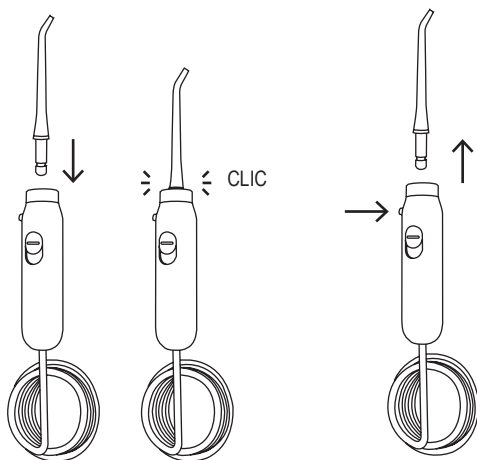
Fig. 2

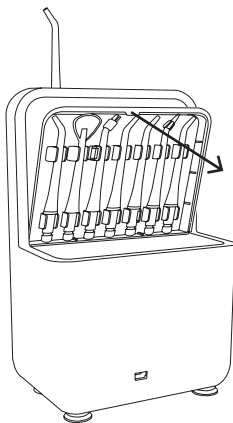**IRRIGADOR BUCAL
HYDRA SMILE****DESCRIPCIÓN**

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English (Original instructions)

**ORAL IRRIGATOR
HYDRA SMILE****DESCRIPTION**

- A 360° nozzle dial
- B Nozzle release button
- C Handle start/stop switch
- D Mode dial (Touch control)
- E ON/OFF selector knob / Pressure regulator
- F Water tank lid
- G Water Tank
- H Nozzle storage panel
- I Supply cord
- J Anti-slip feet
- K Spray Nozzles
 - K1 Standard nozzle (x4)
 - K2 Periodontal nozzle
 - K3 Orthodontic nozzle
 - K4 Tongue cleaner
- L Accessories box

If the model of your appliance does not have the accessories described above, they can also be bought separately from the Technical Assistance Service.

USE AND CARE

- Fully extend the appliance's supply cord before each use.
- Do not use the appliance if the parts or accessories are not properly fitted.

- Do not use the appliance if the accessories attached to it are defective. Replace them immediately.
- Do not use the appliance if it is empty.
- Do not switch the appliance on without water.
- Do not use the appliance if the on/off knob does not work.
- Respect the MAX level (600 ml).
- If you use the appliance in a bathroom or similar place, unplug the appliance from the mains when it is not in use, even if it is only for a short while, as proximity to water involves a risk, even if the appliance is disconnected.
- Disconnect the appliance from the mains when not in use and before undertaking any cleaning task.
- Unplug the appliance from the mains before refilling the water tank.
- Unplug the appliance from the mains before replacing any accessory.
- Keep the appliance in good condition. Check that the moving parts are not misaligned or jammed and make sure there are no broken parts or anomalies that may prevent the appliance from operating correctly.
- Use the appliance and its accessories and tools in accordance with these instructions, taking into account the working conditions and the work to be performed. Use of the appliance for purposes other than those intended could result in a hazardous situation.

- Do not use the appliance with pets or animals.
- Do not direct water under the tongue, into the ear, nose, or other delicate areas.
- Do not use while wearing any oral jewelry, remove prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using the appliance.

INSTALLATION

- Before use, please secure the product to a stable surface using reliable methods such as clamps or screws to prevent it from falling into the water. Drill four holes in the surface with 3mm screws, insert the product on the screws, matching the four holes in the bottom of the appliance.

INSTRUCTIONS FOR USE

BEFORE USE

- Make sure that all the product's packaging has been removed.
- Please read carefully the "Safety advice and warnings" booklet before first use.

FILLING WITH WATER

- You must fill the water tank before turning on the appliance.
- There are two ways of filling the water.

METHOD 1

- Remove the appliance's water tank lid.
- Fill the water tank taking care to observe the MAX level (600 ml).
- Put the water tank lid back into its place.

METHOD 2

- Remove the appliance's water tank from the base.
- Remove the appliance's water tank lid.
- Fill the water tank taking care to observe the MAX level (600 ml).
- Put the water tank lid back into its place.
- Put the water tank back into its housing, ensuring that it is inserted properly.

INSERTING AND REMOVING NOZZLES

- Insert the nozzle into the 360° nozzle dial and press down until you hear a click.
- To remove the nozzle from the handle, press the nozzle release button and pull out the nozzle.

SPRAY NOZZLES

STANDARD NOZZLE (K1)

- Start rinsing from the molars towards the incisors.
- Let the spray head lightly shift along the gums, and pause briefly between teeth.
- Repeat with each side of the teeth.

PERIODONTAL NOZZLE (K2)

- The periodontal nozzle delivers water or anti-bacterial solutions deep into periodontal pockets.
- Select the lowest pressure setting (number 1 and SOFT mode)
- Lean low over your sink and place the soft nozzle against a tooth at a 45-degree angle and gently place the tip under the gum line inserting the nozzle between teeth.
- Continue tracing along the gum line until you finish cleaning.

ORTHODONTIC NOZZLE (K3)

- This nozzle is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work.
- Lean low over the sink and use the brush bristles to clean around dental work, water will flow through the brush to rinse away any debris.

TONGUE CLEANER (K4)

- Lean low over your sink, place the tongue cleaner on your tongue and start the appliance on the lowest pressure setting, increase as comfortable.
- Place the tongue cleaner in the center of your tongue about half way back and pull forward with light pressure.

NOZZLE STORAGE PANEL (FIG.2)

- The nozzle storage panel can be taken out by pulling out the buckle.
- You can take a nozzle out from the panel early before wash and change

the nozzle more conveniently. Or storage the nozzle in the accessory box.

USE

- Extend the supply cord completely before plugging it in.
- Connect the appliance to the mains.
- Guide the appliance to direct the flow of water in the desired direction. The direction of the water flow can be controlled by rotating the the 360° nozzle dial.
- Switch the unit on by selecting the ON setting on the ON/OFF selector knob.
- Select the desired mode in the mode touch control between: NORMAL, SOFT and PULSE.
- To start the water flow, turn on the start/ stop switch.
- When the water flow and the gums are at a 90-degree angle, slightly close you mouth in order to avoid the water spraying out, yet let the water flow from the corner of the mouth into de sink.

STAND-BY FUNCTION

- In order to save energy, the appliance passes to stand-by after 3 minutes.
- To turn it on again, you simply have to switch the ON/OFF selector knob back to the OFF setting, then turn it in a clockwise direction again.

ONCE YOU HAVE FINISHED USING THE APPLIANCE

- Stop the appliance by using the start/ stop switch.

- Switch off the appliance by selecting the OFF position on the ON/OFF selector knob. (E).
- Unplug the appliance from the mains.
- Remove the water from the inside of the appliance.
- Clean the appliance.

CLEANING

- Disconnect the appliance from the mains and allow it to cool before undertaking any cleaning task.
- Mouthwash or antibacterial solutions can be added to this water flosser. After using any special solution, rinse the appliance to prevent clogging by partially filling the reservoir and hose with warm water and running the appliance with the tip pointing into the sink until it is empty.
- Clean the equipment with a damp cloth with a few drops of washing-up liquid and then dry.
- Do not use solvents, or products with an acid or base pH such as bleach, or abrasive products, for cleaning the appliance.
- Never immerse the appliance in water or any other liquid or place it under running water.
- Never submerge the connector in water or any other liquid or place it under running water.

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